

AD-A105 978

LETTERMAN ARMY INST OF RESEARCH PRESIDIO OF SAN FRANC--ETC F/G 6/20
THE MUTAGENIC POTENTIAL OF 4 NITROPHENYL BIS(2-THIENYL)-PHOSPHI--ETC(U)
SEP 81 L J SAUERS, F R PULLIAM, J T FRUYN

UNCLASSIFIED

LAIR-104

NL

For I
40 3
108978

END
DATE
FILMED
11-81
DTIC

LEVEL

13

AD A105978

INSTITUTE REPORT NO. 104

THE MUTAGENIC POTENTIAL OF: 4-nitrophenyl bis(2-thienyl) phosphinate
4-nitrophenyl 2-furyl(methyl) phosphinate
4-cyanophenyl bis(2-furyl) phosphinate
4-nitrophenyl bis(2-furyl) phosphinate

LEONARD J. SAUERS, BA, SP5
FREDDICA R. PULLIAM, BS, SSG
and
JOHN T. FRUIN, DVM, PhD, LTC VC


TOXICOLOGY GROUP,
DIVISION OF RESEARCH SUPPORT

DTIC
ELECT
S OCT 22 1981
A

DTIC FILE COPY

SEPTEMBER 1981

Toxicology Series 17

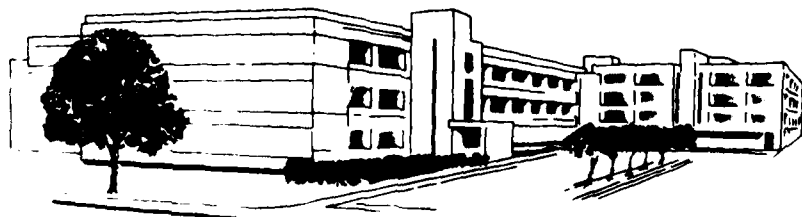

LETTERMAN ARMY INSTITUTE OF RESEARCH PRESIDIO OF SAN FRANCISCO CALIFORNIA 94129

85 10 21

Toxicology Series 17

Reproduction of this document in whole or in part is prohibited except with the permission of the Commander, Letterman Army Institute of Research, Presidio of San Francisco, California 94129. However, the Defense Technical Information Center is authorized to reproduce the document for United States Government purposes.

Destroy this report when it is no longer needed. Do not return it to the originator.

Citation of trade names in this report does not constitute an official endorsement or approval of the use of such items.

This material has been reviewed by Letterman Army Institute of Research and there is no objection to its presentation and/or publication. The opinions or assertions contained herein are the private views of the author(s) and are not to be construed as official or as reflecting the views of the Department of the Army or the Department of Defense. (AR 360-5)

James M. M... J. S. Searcy
(Signature and date)

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

| REPORT DOCUMENTATION PAGE | | READ INSTRUCTIONS BEFORE COMPLETING FORM |
|---|-----------------------|--|
| 1. REPORT NUMBER | 2. GOVT ACCESSION NO. | 3. RECIPIENT'S CATALOG NUMBER |
| LAIR Institute Report No. 104 | AD-4705977 | |
| 4. TITLE (and Subtitle) The Mutagenic Potential of 4 nitro-phenyl bis(2-thienyl)phosphinate; 4-nitrophenyl 2-furyl(methyl)phosphinate; 4-cyanophenyl bis(2-furyl)phosphinate; 4-nitrophenyl bis(2-furyl)phosphinate. | | 5. TYPE OF REPORT & PERIOD COVERED Final Rept. June - September 1981 |
| 7. AUTHOR(s) Léonard J. Sauers, BA, SP5; Freddica R. Pulliam, BS, SSG; John T. Fruin, DVM, PhD, LTC VC; | | 6. PERFORMING ORG. REPORT NUMBER |
| 9. PERFORMING ORGANIZATION NAME AND ADDRESS Toxicology Group, Div. of Research Support Letterman Army Institute of Research Presidio of San Francisco, CA 94129 | | 8. CONTRACT OR GRANT NUMBER(s) (16) 35162772, A 875 |
| 11. CONTROLLING OFFICE NAME AND ADDRESS U.S. Army Medical Research and Development Command Fort Detrick Frederick, MD 21701 | | 10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS Project 3516772A875 WU 304 |
| 14. MONITORING AGENCY NAME & ADDRESS (If different from Controlling Office) (14) LAIR-104 | | 12. REPORT DATE September 1981 |
| | | 13. NUMBER OF PAGES 38 (12) 37 |
| | | 15. SECURITY CLASS. (of this report) UNCLASSIFIED |
| | | 15a. DECLASSIFICATION/DOWNGRADING SCHEDULE |
| 16. DISTRIBUTION STATEMENT (of this Report) APPROVED FOR PUBLIC RELEASE: DISTRIBUTION UNLIMITED. | | |
| 17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report) | | |
| 18. SUPPLEMENTARY NOTES | | |
| 19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Mutagenicity; Toxicology, Ames Assay, 4-nitrophenyl bis(2-thienyl)phosphinate; 4-nitrophenyl 2-furyl(methyl)phosphinate; 4-cyanophenyl bis(2-furyl)phosphinate; 4-nitrophenyl bis(2-furyl)phosphinate. | | |
| 20. ABSTRACT (Continue on reverse side if necessary and identify by block number) The mutagenic potential of 4 nitrophenyl bis(2-thienyl)phosphinate (41*); 4-nitrophenyl 2-furyl(methyl)phosphinate (72*) 4 cyanophenyl bis(2-furyl)phosphinate (82*); 4-nitrophenyl bis(2-furyl)phosphinate (87*) was assessed by using the Ames Salmonella/Mammalian Microsome Mutagenicity Assay. Tester strains TA 98, TA 100, TA 1535, TA 1537 and TA 1538 were exposed to doses ranging from 1 mg/plate to 3.2×10^{-4} mg/plate. It was determined that none of the tested substances had mutagenic potential. *Code number of compound. | | |

DD FORM 1 JAN 73 1473

EDITION OF 1 NOV 65 IS OBSOLETE

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

00032

ABSTRACT

The mutagenic potential of 4 nitrophenyl bis(2-thienyl)phosphinate (41*); 4-nitrophenyl 2-furyl(methyl)phosphinate (72*); 4-cyanophenyl bis(2-furyl)phosphinate (82*); 4-nitrophenyl bis(2-furyl)phosphinate (87*) was assessed by using the Ames Salmonella/Mammalian Microsome Mutagenicity Assay. Tester strains TA 98, TA 100, TA 1535, TA 1537, and 1538 were exposed to doses ranging from 1 mg/plate to 3.2×10^{-4} mg/plate. It was determined that none of the tested substances had mutagenic potential.

* Code number for compound.

| | |
|--------------------|-------------------------------------|
| Accession For | |
| NTIS GRA&I | <input checked="" type="checkbox"/> |
| DTIC TAB | <input type="checkbox"/> |
| Unannounced | <input type="checkbox"/> |
| Justification | |
| By | |
| Distribution/ | |
| Availability Codes | |
| Dist | Avail and/or Special |
| A | |

PREFACE

AMES ASSAY REPORT:

| SUBSTANCE | CODE NO. |
|--|----------|
| 4 nitrophenyl bis(2-thienyl)phosphinate | 41 |
| 4-nitrophenyl 2-furyl(methyl)phosphinate | 72 |
| 4-cyanophenyl bis(2-furyl)phosphinate | 82 |
| 4-nitrophenyl bis(2-furyl)phosphinate | 87 |

TESTING FACILITY: Letterman Army Institute of Research
Presidio of San Francisco, CA 94129

SPONSOR: Biomedical Laboratory, Aberdeen Proving Grounds
Aberdeen, MD 21005

PROJECT: Toxicity Testing of Phosphinate Compounds - 35162772A875

GLP STUDY NUMBER: 81013

STUDY DIRECTOR: LTC John T. Fruin D.V.M., PhD.

CO-PRINCIPAL INVESTIGATORS: SSG Freddica R. Pulliam, B.S.
SP5 Leonard J. Sauers, B.A.

RAW DATA: A copy of the final report, study protocol and retired SOPs will be maintained in the LAIR archives. Test substances were provided by sponsor. Chemical, analytical, stability, purity, etc. data are available from the sponsor.


PURPOSE: To determine the mutagenic potential of the above compounds using the Ames Assay. Tester strains TA 98, TA 100, TA 1535, TA 1537, and TA 1538 were used.


ACKNOWLEDGMENTS


The authors wish to thank John Dacey and SP4 Larry Mullen, BS for their assistance in performing the research and for help in preparation of this report.

Signatures of Principal Scientists
Involved in the Study

We, the undersigned, believe the study, GLP number 81013, described in this report to be scientifically sound and the results and interpretation to be valid. The study was conducted to comply to the best of our ability with the Good Laboratory Practice Regulations outlined by the Food and Drug Administration.

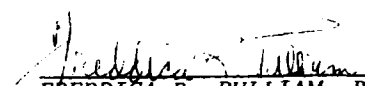
FREDDICA R. PULLIAM, BS Date 8 Jul 81
SSG
Co-Investigator


JOHN T. FRUIN, DVM, PhD Date 7 July 81
LTC, VC
Study Director


LEONARD J. SAUERS, BA Date 9 July 81
SP5
Co-Investigator



DEPARTMENT OF THE ARMY
LETTERMAN ARMY INSTITUTE OF RESEARCH
PRESIDIO OF SAN FRANCISCO, CALIFORNIA 94129

REPLY TO
ATTENTION OF:

NRD-UL-4-7A

21 July 1981

MEMORANDUM FOR RECORD

SUBJECT: Report of GLP Compliance

I hereby certify that in relation to LAIR GLP study 81013 the following inspection was made:

17 June 1981

Routine inspections with no adverse findings are reported quarterly, thus this inspection is also included in the July 1981 report to management.

JOHN C. JOHNSON
CPT, MS
Quality Assurance Officer

TABLE OF CONTENTS

| | |
|---|-----|
| Abstract..... | i |
| Preface..... | iii |
| Acknowledgments..... | iv |
| Signatures of Principal Scientists..... | v |
| Report of Quality Assurance Unit..... | vi |
| Table of Contents..... | vii |

BODY OF REPORT

INTRODUCTION

| | |
|--|---|
| Rationale for using the Ames Assay..... | 1 |
| Description of Test, Rationale for strain selection..... | 1 |
| Description of Strains, History, Methods, and Data..... | 2 |

METHODS

| | |
|--|---|
| Rationale for Dosage Levels and Response Tabulations.... | 3 |
| Test Format..... | 3 |
| Statistical Analysis..... | 4 |

| | |
|-----------------------------|---|
| RESULTS AND DISCUSSION..... | 4 |
|-----------------------------|---|

| | |
|-----------------|---|
| CONCLUSION..... | 5 |
|-----------------|---|

| | |
|---------------------|---|
| RECOMMENDATION..... | 5 |
|---------------------|---|

| | |
|-----------------|---|
| REFERENCES..... | 6 |
|-----------------|---|

| | |
|------------------------------------|---|
| APPENDIX (Tables 1 through 6)..... | 7 |
|------------------------------------|---|

| | |
|------------------------|----|
| DISTRIBUTION LIST..... | 28 |
|------------------------|----|

Rationale for using the Ames Assay

The Ames Salmonella/Mammalian Microsome Mutagenicity Test is one of a standard bank of tests used by our laboratory for the assessment of the mutagenic potential of a test substance. It is a short-term screening assay for the prediction of potential mutagenic agents in mammals. It is inexpensive when compared to in vivo tests, yet is highly predictive and reliable in its ability to detect mutagenic activity and therefore carcinogenic probability (1). It relies on basic genetic principles and allows for the incorporation of a mammalian microsome enzyme system to increase sensitivity through enzymatically altering the test substance into an active metabolite. It has proven highly effective in assessing human risk (1).

Description of Test (Rationale for the selection of strains)

The test was developed by Bruce Ames, Ph.D. from the University of California-Berkeley. The test involves the use of several different genetically altered strains of Salmonella typhimurium, each with a specific mutation in the histidine operon (2). The test substance demonstrates mutagenic potential if it is able to revert the mutation in the bacterial histidine operon back to the wild type and thus reestablish prototrophic growth within the test strain. This reversion also can occur spontaneously due to a random mutational event. If, after adding a test substance, the number of revertants is significantly greater than the spontaneous reversion rate, then the test substance physically altered the locus involved in the operon's mutation and is able to induce point mutations and genetic damage (2).

In order to increase the sensitivity of the test system, two other mutations in the Salmonella are used (2). To insure a higher probability of uptake of test substance, the genome for the lipopolysacchride layer (LP) is mutated and allows larger molecules to enter the bacteria. Each strain has another induced mutation which causes loss of excision repair mechanisms. Since many chemicals are not by themselves mutagenic but have to be activated by an enzymatic process, a mammalian microsome system is incorporated. These microsomal enzymes are obtained from livers of rats induced with Aroclor 1254; the enzymes allow for the expression of the metabolites in the mammalian system. This activated rat liver microsomal enzyme homogenate is termed S-9.

Description of Strains (History of the strains used, methods to monitor the integrity of the organisms, and data pertaining to current and historical controls and spontaneous reversion rates)

The test consists of using five different strains of Salmonella typhimurium that are unable to grow in absence of histidine because of a specific mutation in the histidine operon. This histidine requirement is verified by attempting to grow the tester strains on minimal glucose agar (MGA) plates, both with and without histidine. The dependence on this amino acid is shown when growth occurs only in its presence. The plasmids in strains TA 98 and TA 100 contain an ampicillin resistant R factor. Strains deficient in this plasmid demonstrate a zone of growth inhibition around an ampicillin impregnated disc. The alteration of the LP layer allows uptake by the Salmonella of larger molecules. If a crystal violet impregnated disc is placed onto a plate containing any one of the bacterial strains, a zone of growth inhibition will occur because the LP layer is altered. The absence of excision repair mechanisms can be determined by using ultraviolet (UV) light. These mechanisms function primarily by repairing photodimers between pyrimidine bases; exposure of bacteria to UV light will activate the formation of these dimers and cause cell lethality, since excision of these photodimers can not be made. The genetic mutation resulting in UV sensitivity also induces a dependence by the Salmonella to biotin. Therefore, this vitamin must be added. In order to prove that the bacteria are responsive to the mutation process, positive controls are run with known mutagens. If after exposure to the positive control substance, a larger number of revertants are obtained, then the bacteria are adequately responsive. Sterility controls are performed to determine the presence of contamination. Sterility of the test compound is also confirmed in each first dilution. Verification of the tester strains occurs spontaneously with the running of each assay. The value of the spontaneous reversion rate is obtained using the same inoculum of bacteria that is used in the assay (3).

Strains were obtained directly from Dr. Ames, University of California, Berkeley, propagated and then maintained at -80 C in our laboratory. Before any substance was tested, quality controls were run on the bacterial strains to establish the validity of their special features and also to determine the spontaneous reversion rate (2). Records are maintained of all the data, to determine if deviations from the set trends have occurred.

We compared the spontaneous reversion values with our own historical values and those cited by Ames et al (2). Our conclusions are based on the spontaneous reversion rate compared to the experimentally induced rate of mutation. When operating effectively, these strains detect substances that cause base pair

mutations (TA 1535, TA 100) and frameshift mutations (TA 1537, TA 1538 and TA 98) (2).

METHODS (3)

Rationale for Dosage Levels and Dose Response Tabulations

To insure readable and reliable results, a sublethal concentration of the test substance had to be determined. This toxicity level was found by using MGA plates, various concentrations of the substance, and approximately 10^8 cells of TA 100 per plate, unless otherwise specified. Top agar containing trace amounts of histidine and biotin were placed on MGA plates. TA 100 is used because it is the most sensitive strain. Strain verification was confirmed on the bacteria, along with a determination of the spontaneous reversion rate. After incubation, the growth was observed on the plates. (The auxotrophic Salmonella will replicate a few times and potentially express a mutation. When the histidine and biotin supplies are exhausted, only those bacteria that reverted to the prototrophic phenotype will continue to reproduce and form macrocolonies; the remainder of the bacteria comprises the background lawn. The minimum toxic level is defined as the lowest serial dilution at which decreased macrocolony formation, below that of the spontaneous revertant rate, and an observable reduction in the density of the background lawn occurs.) A maximum dose of 1 mg/plate is used when no toxicity is observed. The densities were recorded as normal slight, and no growth.

Test Format

After we validated our bacterial strains and determined the optimal dosage of the test substance, we began the Ames Assay. In the actual experiment, 0.1ml of the particular strain of Salmonella (10^8 cells) and the specific dilutions of the test substance were added to 2 ml of molten top agar, which contained trace amounts of histidine and biotin. Since survival is better from cultures which have just passed the log phase, the Salmonella strains were used 16 hours (maximum) after initial inoculation into nutrient broth. The dose of the test substance spanned more than a 1000- fold, decreasing from the minimum toxic level by a dilution factor of 5. All the substances were tested with and without S-9 microsome fraction. The S-9 mixture which was previously titrated at an optimal strength was added to the molten top agar. After all the ingredients were added, the top agar was vortexed, then overlaid on minimum glucose agar plates. These plates contained 2% glucose and Vogel Bonner "E" Concentrate (4). The water used in this medium and all reagents came from a polymetric system. Plates were incubated, upside down in the dark at 37 C for 48 hours. Plates were prepared in triplicate and the average revertant counts were recorded. The corresponding number

of revertants obtained was compared to the number of spontaneous revertants; the conclusions were recorded statistically. A correlated dose response is considered necessary to declare a substance as a mutagen. Commoner (5), in his report, "Reliability of Bacterial Mutagenesis Techniques to Distinguish Carcinogenic and Non-Carcinogenic Chemical," and McCann et al (1) in their paper, "Detection of Carcinogens as Mutagen: Assay of over 300 Chemicals," have concurred on the test's ability to detect mutagenic potential.

Statistical Analysis

Quantitative evaluation was ascertained by two independent methods. Ames et al (2) assumed that a compound which caused twice the spontaneous reversion rate is mutagenic. Commoner (5), developed the MUTAR Ratio, which is stated in the following equation:

$$\text{MUTAR} = (E - C)/C_{AV}$$

Here, C is the number of spontaneous revertant colonies on control plates obtained on the same day and with the same treatment and strains. E is the number of revertants in response to the compound; C_{AV} is the number of spontaneous revertants on control plates calculated from historical records. The explanation of the results of this equation can be determined by the method of Commoner (5). This variation determines the probability of correctly classifying substances as carcinogens on the basis of their mutagenic activity. The E values were recorded by strain, with and without S-9. Values for C and C_{AV} were recorded separately.

We used the formula and logged all values for our permanent records.

RESULTS AND DISCUSSION

Throughout this report, each of the test substances will be referred to by the respective code number:

| <u>Substance</u> | <u>Code No.</u> |
|--|-----------------|
| 4 nitrophenyl bis(2-thienyl)phosphinate | 41 |
| 4-nitrophenyl 2-furyl(methyl)phosphinate | 72 |
| 4-cyanophenyl bis(2-furyl)phosphinate | 82 |
| 4-nitrophenyl bis(2-furyl)phosphinate | 87 |

On 1 June 1981, the Toxicity Level Determination was performed on the 4 test chemicals. All positive, negative and sterility controls for this experiment were normal (Table 1). At the highest dose used, 1.0 mg/plate, no toxicity was observed (Table 2A-2D).

On 17 June 1981, the Ames Assay was performed using the 4 test substances. For this experiment, all sterility and strain verification controls were normal (Table 3). Expected responses were observed for all negative and positive controls (Table 4).

For all the chemicals tested, there were no incidence of mutagenicity (Table 5A-5D).

The MUTAR values listed in Tables 6A-6D were within the normal limits.

CONCLUSION

On the basis of the Ames Assay, test compounds 41, 72, 82, and 87 are not mutagenic at the levels tested.

RECOMMENDATION

We recommend that organophosphinate compounds 41, 72, 82, and 87 be tested by using other toxicological testing systems if efficacy tests show these chemicals to be promising antidotes.

REFERENCES

1. McCANN, J., E. CHOI, E. YAMASAKI, and B. N. AMES. Detection of carcinogens as mutagens in the Salmonella/microsome test: Assay of 300 chemicals. Proc Nat Acad Sci, USA 72:5135-5139, 1975
2. AMES, B. N., J. McCANN and E. YAMASAKI. Methods for detection of carcinogens and mutagens with Salmonella/mammalian microsome mutagenicity test. Mutation Res 31: 347-364, 1975
3. LAIR SOP OP-STX-1, Ames Salmonella/mammalian microsome mutagenicity test, 1 March 1981
4. VOGEL, H. J. and D. M. BONNER. Acetylornithinase of E. coli: Partial purification and some properties. J Biol Chem 218: 97-106, 1956
5. COMMONER, B. Reliability of the bacterial mutagenesis techniques to distinguish carcinogenic and non-carcinogenic chemicals. EPA 600/1 76-022, 1976

LIST OF TABLES

| | Code | Page |
|--|------|------|
| Table 1 Strain Verification for Toxicity Level Determination | | 9 |
| Table 2A Toxicity Level Determination | 41 | 10 |
| Table 2B Toxicity Level Determination | 72 | 11 |
| Table 2C Toxicity Level Determination | 82 | 12 |
| Table 2D Toxicity Level Determination | 87 | 13 |
| Table 3 Strain Verification and Sterility Controls | | 14 |
| Table 4 Positive and Negative Controls | | 15 |
| Table 5A Ames Assay Worksheet | 41 | 16 |
| Table 5B Ames Assay Worksheet | 72 | 18 |
| Table 5C Ames Assay Worksheet | 82 | 20 |
| Table 5D Ames Assay Worksheet | 87 | 22 |
| Table 6A Mutagenic Activity Ratio | 41 | 24 |
| Table 6B Mutagenic Activity Ratio | 72 | 25 |
| Table 6C Mutagenic Activity Ratio | 82 | 26 |
| Table 6D Mutagenic Activity Ratio | 87 | 27 |

APPENDIX

TABLE 1
STRAIN VERIFICATION FOR TOXICITY LEVEL DETERMINATION
Salmonella/Microsome Assay

| Strain No. | Histidine Requirements | Ampicillin Resistance | uvr-B Deletion | rfa Crystal Violet | Sterility Control | Response (a) |
|--|------------------------|-----------------------|----------------|--------------------|-------------------|--------------|
| TA 100 | NG | G | NG | 15.46 mm | NG | + |
| TA 1537 | NG | NG | NG | 14.11 mm | NG | + |
| WT | G | NA | G | NA | NA | + |
| Diluent | NA | NA | NA | NA | NG | + |
| Positive Control - MNNG - Average - 1612 | | | | | | |
| Test Compound (s) | | | | | | |
| (a) 41 | NA | NA | NA | NA | NG | + |
| (b) 72 | NA | NA | NA | NA | NG | + |
| (c) 82 | NA | NA | NA | NA | NG | + |
| (d) 87 | NA | NA | NA | NA | NG | + |
| (e) NA | NA | NA | NA | NA | NA | NA |

G = Growth; NG = No Growth; NT = Not Tested; NA = Not Applicable;
WT = Wild Type; (a) + = Expected Response; - = Unexpected Response

Spontaneous Revertants

| Strain | Time | | | | Average |
|--------|-----------|-----|-----|-----|---------|
| TA 100 | Beginning | 146 | 148 | 155 | 140 |
| TA 100 | End | 122 | 118 | 149 | |

Test Inoculated By: Sauers, Pulliam, Dacey, Mullen Date 1 June 1981

Test Read By: Sauers, Pulliam Date 3 June 1981

TABLE 2A
TOXICITY LEVEL DETERMINATION
Salmonella/Microsome Assay

Substance assayed: (1) Code #41 (2) _____

(3) _____ (4) _____ (5) _____

Date: 3 June 1981 Performed by: Sauers, Pulliam, Dacey, Mullen

Substance dissolved in: (1) DMSO (2) _____ (3) _____

(4) _____ (5) _____

Visual estimation of background lawn on
Nutrient Agar Plates: NG = no growth
ST = slight growth
NL = normal growth

TA 100
Revertant Plate Count

| Test Compound Concentration | Plate #1 | Plate #2 | Plate #3 | Average | Background Lawn |
|--------------------------------|----------|----------|----------|---------|--------------------|
| 1.0 mg/pl | 135 | 105 | 138 | 126 | NL |
| 10 ⁻¹ | 125 | 150 | 106 | 127 | NL |
| 10 ⁻² | 121 | 124 | 117 | 121 | NL |
| 10 ⁻³ | 138 | 132 | 126 | 132 | NL |
| 10 ⁻⁴ | 128 | 118 | 132 | 126 | NL |
| 10 ⁻⁵ | 134 | 125 | 152 | 137 | NL |
| 10 ⁻⁶ | 154 | 139 | 123 | 139 | NL |
| 10 ⁻⁷ | 139 | 160 | 179 | 159 | NL |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |

TABLE 2B
TOXICITY LEVEL DETERMINATION
Salmonella/Microsome Assay

Substance assayed: (1) Code #72 (2) _____

(3) _____ (4) _____ (5) _____

Date: 3 June 1981 Performed by: Pulliam, Sauers, Dacey, Mullen

Substance dissolved in: (1) DMSO (2) _____ (3) _____

(4) _____ (5) _____

Visual estimation of background lawn on
Nutrient Agar Plates: NG = no growth
ST = slight growth
NL = normal growth

TA 100
Revertant Plate Count

| Test Compound Concentration | Plate #1 | Plate #2 | Plate #3 | Average | Background Lawn |
|--------------------------------|----------|----------|----------|---------|--------------------|
| 1.0 mg/pl | 180 | 194 | 166 | 180 | NL |
| 10 ⁻¹ | 193 | 232 | 228 | 218 | NL |
| 10 ⁻² | 206 | 193 | 181 | 193 | NL |
| 10 ⁻³ | 194 | 193 | 190 | 192 | NL |
| 10 ⁻⁴ | 130 | 135 | 113 | 126 | NL |
| 10 ⁻⁵ | 125 | 116 | 89 | 110 | NL |
| 10 ⁻⁶ | 108 | 113 | 142 | 121 | NL |
| 10 ⁻⁷ | 159 | 140 | 146 | 148 | NL |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |

TABLE 2C
TOXICITY LEVEL DETERMINATION
Salmonella/Microsome Assay

Substance assayed: (1) Code #82 (2) _____

(3) _____ (4) _____ (5) _____

Date: 3 June 1981 Performed by: Sauers, Pulliam, Dacey, Mullen

Substance dissolved in: (1) DMSO (2) _____ (3) _____

(4) _____ (5) _____

Visual estimation of background lawn on
Nutrient Agar Plates: NG = no growth
ST = slight growth
NL = normal growth

TA 100
Revertant Plate Count

| Test Compound Concentration | Plate #1 | Plate #2 | Plate #3 | Average | Background Lawn |
|--------------------------------|----------|----------|----------|---------|--------------------|
| 1.0 mg/pl | 141 | 148 | 108 | 132 | NL |
| 10 ⁻¹ | 108 | 134 | 154 | 132 | NL |
| 10 ⁻² | 132 | 102 | 142 | 125 | NL |
| 10 ⁻³ | 133 | 148 | 144 | 142 | NL |
| 10 ⁻⁴ | 130 | 143 | 164 | 146 | NL |
| 10 ⁻⁵ | 154 | 125 | 123 | 134 | NL |
| 10 ⁻⁶ | 136 | 123 | 109 | 123 | NL |
| 10 ⁻⁷ | 139 | 132 | 139 | 137 | NL |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |

TABLE 2D
TOXICITY LEVEL DETERMINATION
Salmonella/Microsome Assay

Substance assayed: (1) Code #87 (2) _____

(3) _____ (4) _____ (5) _____

Date: 3 June 1981 Performed by: Sauers, Pulliam, Dacey, Mullen

Substance dissolved in: (1) DMSO (2) _____ (3) _____

(4) _____ (5) _____

Visual estimation of background lawn on
Nutrient Agar Plates: NG = no growth
ST = slight growth
NL = normal growth

TA 100
Revertant Plate Count

| Test Compound Concentration | Plate #1 | Plate #2 | Plate #3 | Average | Background Lawn |
|--------------------------------|----------|----------|----------|---------|--------------------|
| 1.0 mg/pl | 125 | 103 | 121 | 116 | NL |
| 10 ⁻¹ | 167 | 158 | 128 | 151 | NL |
| 10 ⁻² | 140 | 123 | 132 | 132 | NL |
| 10 ⁻³ | 140 | 158 | 131 | 143 | NL |
| 10 ⁻⁴ | 141 | 119 | 136 | 132 | NL |
| 10 ⁻⁵ | 139 | 97 | 123 | 120 | NL |
| 10 ⁻⁶ | 139 | 165 | 116 | 140 | NL |
| 10 ⁻⁷ | 122 | 140 | 182 | 148 | NL |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |
| | | | | | |

TABLE 3
STRAIN VERIFICATION CONTROL

| Strains | Histidine Requirement | Ampicillin Resistance | UV | Sensitivity to Crystal Violet | Sterility Control | Response (1) |
|---------|-----------------------|-----------------------|----|-------------------------------|-------------------|--------------|
| 98 | NG | G | NG | 14.15 mm | NG | + |
| 100 | NG | G | NG | 14.51 mm | NG | + |
| 1535 | NG | NA | NG | 13.89 mm | NG | + |
| 1537 | NG | NG | NG | 14.75 mm | NG | + |
| 1538 | NG | NA | NG | 15.11 mm | NG | + |
| WT | G | NA | G | NA | NA | + |

STERILITY CONTROL

His-Bio Mix Initial: NG End: NG Diluent: NG
 Top Agar Initial: NG End: NG MGA Plate: NG
 S-9 Mix Initial: NG End: NG Nutrient Broth: NG
 Test Compound (a) 41-NG (b) 72-NG (c) 82-NG (d) 87-NG (e) NA (f) NA

G = Growth NG = No Growth NT = Not Tested NA = Not Applicable WT = Wild Type

Study Number: 81013 By: Sauers, Pulliam,

Date: 17 June 1981 Dacey, Mullen (1) + = expected response
 - = unexpected response

TABLE 4

| <u>SPONTANEOUS REVERTANT RATE AND POSITIVE CONTROL REVERTANT RATE</u> | | | | | |
|---|-----------------------------------|----------------------|----------------------------------|-------------------------------------|---|
| <u>Compd.</u> | <u>Amount of Compd. Added</u> | <u>S-9 Added</u> | <u>98</u> | <u>Strain Number</u> | |
| | | | | <u>1535</u> | <u>1537</u> |
| | | | | <u>100</u> | <u>1538</u> |
| AF | 2 ug/plate | yes | (372,593,430) (465) | (347,366,471) (395) | (212,452,469) (378) |
| BF | 2 ug/plate | yes | (88,91,153) (111) | (186,207,306) (233) | (63,28,60) (50) |
| DMBA | 20 ug/plate | yes | (70,39,50) (53) | (266,207,249) (241) | (36,37,13) (29) |
| MNNG | 2 ug/plate | no | | (369,312,425) (369) | |
| | 20 ug/plate | no | | (87,102,NG) (95) | |
| <u>Strain Performance</u> | | | | | |
| Spontaneous Revertants | | | | | |
| | before after | no | (23,15,17) (18,18,18) (17) | (99,86,82) (85,69,92) (86) | (7,12,8) (15,21,9) (12) |
| | | | | | (4,7,7) (6,13,7) (7) |
| | before after | yes | (18,17,20) (15,12,11) (16) | (87,101,105) (75,109,85) (94) | (10,12,9) (17,15,10) (13) |
| | | | | | (4,4,6) (4,7,9) (6) |
| | | | | | (9,9,10) (15,15,18) (13) (6,7,9) (20,11,10) (11) |

Study Number: 81013

Date: 17 Jun 31 By: Savers, Pulliam, Dacey, Mullen

TABLE 5A

NUMBER OF REVERTANTS/PLATE

| Compd. | Amount of Compd. Added | S-9 Added | Strain Number | | | |
|----------|---------------------------|--------------|--------------------|------------------------|--------------------------------|--------------------------------|
| | | | 98 | 100 | 1535 | 1537 |
| Code #41 | 1 mg/plate | no | (Toxic) (Toxic) | (Toxic) (Toxic) | (Toxic,Toxic,Toxic) (Toxic) | (Toxic,Toxic,Toxic) (Toxic) |
| | | yes | (18,32,21) (24) | (125,97,140) (121) | (9,10,9) (9) | (7,5,4) (5) |
| Code #41 | 0.2 mg/plate | no | (18,15,10) (14) | (85,100,88) (91) | (14,12,9) (12) | (3,4,5) (4) |
| | | yes | (18,28,27) (24) | (123,144,129) (132) | (10,12,12) (11) | (3,6,7) (5) |
| Code #41 | 0.04 mg/plate | no | (14,14,14) (14) | (98,71,95) (88) | (12,13,16) (14) | (6,8,5) (6) |
| | | yes | (18,19,28) (22) | (117,110,110) (112) | (7,9,10) (9) | (6,6,10) (7) |

-continued

TABLE 5A, concluded

| NUMBER OF REVERTANTS/PLATE | | | | | | |
|----------------------------|---------------------------|--------------|--------------------|------------------------|--------------------|-----------------|
| Compd. | Amount of Compd. Added | S-9 Added | 98 | 100 | Strain Number | |
| | | | | | 1535 | 1537 |
| Code #41 | 0.008 mg/plate | no | (20,23,24) (22) | (117,126,135) (126) | (17,13,17) (16) | (3,8,8) (6) |
| | | yes | (20,14,11) (15) | (93,90,113) (99) | (12,10,8) (18) | (6,5,10) (7) |
| Code #41 | 0.0016 mg/plate | no | (11,11,21) (14) | (97,84,93) (91) | (12,11,6) (10) | (4,9,5) (6) |
| | | yes | (29,17,17) (21) | (133,120,126) (126) | (11,7,10) (9) | (7,6,6) (6) |
| Code #41 | 0.00032 mg/plate | no | (12,12,21) (15) | (98,101,87) (95) | (17,20,14) (17) | (5,9,9) (8) |
| | | yes | (25,16,37) (26) | (105,80,119) (101) | (20,18,10) (16) | (11,8,7) (9) |

Study Number: 81013

Date: 17 June 81

By: Sauer, Pulliam, Dacey, Mullen

TABLE 5B

NUMBER OF REVERTANTS/PLATE

| Compd. | Amount of Compd. Added | S-9 Added | Strain Number | | | |
|----------|---------------------------|--------------|--------------------|------------------------|--------------------|-----------------|
| | | | 98 | 100 | 1535 | 1537 |
| Code #72 | 1 mg/plate | no | (7,9,15) (10) | (82,91,72) (82) | (6,8,14) (9) | (6,9,7) (7) |
| | | yes | (23,17,18) (19) | (107,111,108) (109) | (10,7,8) (8) | (6,7,14) (9) |
| Code #72 | 0.2 mg/plate | no | (13,14,18) (15) | (68,62,65) (65) | (11,12,12) (12) | (4,5,5) (5) |
| | | yes | (19,18,12) (16) | (124,138,110) (124) | (8,12,5) (8) | (8,6,5) (6) |
| Code #72 | 0.04 mg/plate | no | (19,11,14) (15) | (79,86,71) (80) | (17,11,17) (15) | (7,5,6) (6) |
| | | yes | (25,24,26) (25) | (101,108,101) (103) | (12,7,12) (10) | (7,5,8) (7) |

-continued

Study Number: 81013

Date: 17 Jun 81

By: Savers, Pulliam, Dacey, Mullen

NUMBER OF REVERTANTS/FLATE

By: Sauers, Pulliam, Dacey, Mullen

TABLE 5C
NUMBER OF REVERTANTS/PLATE

| Compd. | Amount of Compd. Added | S-9 Added | 98 | 100 | Strain Number | | |
|----------|---------------------------|--------------|--------------------|------------------------|--------------------|-----------------|--------------------|
| | | | | | 1535 | 1537 | 1538 |
| Code #82 | 1 mg/plate | no | (10,14,10) (11) | (74,70,70) (71) | (10,11,15) (12) | (4,4,7) (5) | (17,10,18) (15) |
| | | yes | (27,16,13) (19) | (91,105,109) (102) | (6,7,11) (8) | (8,4,7) (6) | (14,12,20) (15) |
| Code #82 | 0.2 mg/plate | no | (19,17,10) (15) | (86,77,54) (72) | (12,6,9) (9) | (5,9,7) (7) | (15,10,17) (14) |
| | | yes | (21,30,30) (27) | (125,109,107) (114) | (11,11,9) (10) | (14,7,4) (8) | (16,19,11) (15) |
| Code #82 | 0.04 mg.plate | no | (12,10,12) (11) | (77,64,56) (66) | (14,25,10) (16) | (3,4,7) (5) | (15,7,10) (11) |
| | | yes | (18,14,15) (16) | (101,79,90) (90) | (9,12,8) (10) | (3,7,7) (6) | (17,20,11) (16) |

-continued

Study Number: 81013 Date: 17 Jun 81 By: Sauters, Pulliam, Dacey, Mullen

TABLE 5C, concluded

NUMBER OF REVERTANTS/PLATE

| Compd. | Amount of Compd. Added | S-9 Added | 98 | 100 | Strain Number | |
|----------|---------------------------|--------------|--------------------|------------------------|--------------------|-----------------|
| | | | | | 1535 | 1537 |
| Code #82 | 0.008 mg/plate | no | (11,17,12) (13) | (82,75,64) (74) | (12,15,12) (13) | (9,5,8) (7) |
| | | yes | (23,14,17) (18) | (126,122,108) (119) | (8,6,8) (7) | (10,9,6) (8) |
| Code #82 | 0.0016 mg/plate | no | (10,19,15) (15) | (89,72,71) (77) | (12,12,9) (11) | (4,7,5) (5) |
| | | yes | (25,20,20) (22) | (77,83,113) (91) | (9,3,18) (12) | (6,10,5) (7) |
| Code #82 | 0.00032 mg/plate | no | (14,9,11) (11) | (96,71,111) (93) | (12,17,12) (14) | (5,6,6) (6) |
| | | yes | (30,29,14) (24) | (80,99,94) (91) | (11,6,8) (8) | (4,6,5) (5) |

Study Number: 81013 Date: 17 Jun 81 By: Sauers, Pulliam, Dacey, Mullen

TABLE 5D
NUMBER OF REVERTANTS/PLATE

| Compd. | Amount of Compd. Added | S-9 Added | Strain Number | | | |
|----------|---------------------------|--------------|--------------------|------------------------|--------------------|-----------------|
| | | | 98 | 100 | 1535 | 1537 |
| Code #87 | 1 mg/plate | no | (10,8,11) (10) | (78,64,98) (80) | (9,11,6) (9) | (6,6,5) (6) |
| | | yes | (35,21,19) (25) | (102,100,134) (112) | (14,15,7) (12) | (10,9,5) (8) |
| Code #87 | 0.2 mg/plate | no | (14,10,10) (11) | (73,95,74) (81) | (9,21,10) (13) | (7,6,10) (8) |
| | | yes | (21,21,18) (20) | (135,91,92) (106) | (17,8,7) (11) | (9,7,10) (9) |
| Code #87 | 0.04 mg/plate | no | (8,12,13) (11) | (83,91,77) (84) | (21,13,20) (18) | (5,4,5) (5) |
| | | yes | (20,18,22) (20) | (79,89,134) (101) | (7,12,12) (10) | (6,12,7) (8) |

-continued

Study Number: 81013 Date: 17 Jun 81 By: Sauers, Pulliam, Dacey, Mullen

NUMBER OF REVERTANTS/PLATE

Study Number: 31013 Date: 17 Jun 81
By: Sauers, Pulliam, Dacev, Mullen

TABLE 6A
MUTAGENIC ACTIVITY RATIO

Substance Assayed: Code #41 Dissolved in: DMSO

Study Number: 81013 Date: 15 July 1981 By: Sauers

| Concentration | Strain | MUTAR (act) | MUTAR | Concentration | Strain | MUTAR (act) | MUTAR |
|----------------|---------|----------------|-------|----------------|---------|----------------|-------|
| 1.0 mg/plate | TA 98 | 0.32 | * | 0.008 mg/plate | TA 1535 | * | 0.26 |
| 0.2 mg/plate | TA 98 | 0.32 | * | 0.0016 mg/pl. | TA 1535 | * | * |
| 0.04 mg/plate | TA 98 | 0.24 | * | 0.00032 mg/pl. | TA 1535 | 0.27 | 0.32 |
| 0.008 mg/plate | TA 98 | * | 0.24 | | | | |
| 0.0016 mg/pl. | TA 98 | 0.2 | * | 1.0 mg/plate | TA 1537 | * | * |
| 0.00032 mg/pl. | TA 98 | 0.4 | * | 0.2 mg/plate | TA 1537 | * | * |
| | | | | 0.04 mg/plate | TA 1537 | 0.15 | * |
| 1.0 mg/plate | TA 100 | 0.25 | * | 0.008 mg/plate | TA 1537 | 0.15 | * |
| 0.2 mg/plate | TA 100 | 0.35 | 0.05 | 0.0016 mg/pl. | TA 1537 | * | * |
| 0.04 mg/plate | TA 100 | 0.17 | 0.02 | 0.00032 mg/pl. | TA 1537 | 0.46 | 0.15 |
| 0.008 mg/plate | TA 100 | 0.05 | 0.42 | | | | |
| 0.0016 mg/pl. | TA 100 | 0.29 | 0.05 | 1.0 mg/plate | TA 1538 | 0.43 | * |
| 0.00032 mg/pl. | TA 100 | 0.06 | 0.1 | 0.2 mg/plate | TA 1538 | 0.53 | * |
| | | | | 0.04 mg/plate | TA 1538 | 0.27 | * |
| 1.0 mg/plate | TA 1535 | * | * | 0.008 mg/plate | TA 1538 | * | 0.21 |
| 0.2 mg/plate | TA 1535 | * | * | 0.0016 mg/pl. | TA 1538 | 0.48 | * |
| 0.04 mg/plate | TA 1535 | * | 0.13 | 0.00032 mg/pl. | TA 1538 | 0.37 | 0.07 |

(act): S-9 fraction was added

* : calculated value resulted in a negative MUTAR or zero MUTAR

TABLE 6B
MUTAGENIC ACTIVITY RATIO

Substance Assayed: Code #72 Dissolved in: DMSO
Study Number: 31013 Date: 15 July 1981 By: Sauers

| Concentration | Strain | MUTAR (act) | MUTAR | Concentration | Strain | MUTAR (act) | MUTAR |
|----------------|---------|----------------|-------|----------------|---------|----------------|-------|
| 1.0 mg/plate | TA 98 | 0.12 | * | 0.008 mg/plate | TA 1535 | * | * |
| 0.2 mg/plate | TA 98 | * | * | 0.0016 mg/pl. | TA 1535 | * | 0.06 |
| 0.04 mg/plate | TA 98 | 9.36 | * | 0.00032 mg/pl. | TA 1535 | * | 0.26 |
| 0.008 mg/plate | TA 98 | * | * | | | | |
| 0.0016 mg/pl. | TA 98 | 0.24 | * | 1.0 mg/plate | TA 1537 | 0.46 | * |
| 0.00032 mg/pl. | TA 98 | 0.16 | * | 0.2 mg/plate | TA 1537 | * | * |
| | | | | 0.04 mg/plate | TA 1537 | 0.15 | * |
| 1.0 mg/plate | TA 100 | 0.14 | * | 0.008 mg/plate | TA 1537 | 0.15 | * |
| 0.2 mg/plate | TA 100 | 0.28 | * | 0.0016 mg/pl. | TA 1537 | * | * |
| 0.04 mg/plate | TA 100 | 0.08 | * | 0.00032 mg/pl. | TA 1537 | * | * |
| 0.008 mg/plate | TA 100 | 0.17 | * | | | | |
| 0.0016 mg/pl. | TA 100 | * | * | 1.0 mg/plate | TA 1538 | * | * |
| 0.00032 mg/pl. | TA 100 | 0.14 | 0.01 | 0.2 mg/plate | TA 1538 | 0.43 | * |
| | | | | 0.04 mg/plate | TA 1538 | 0.48 | * |
| 1.0 mg/plate | TA 1535 | * | * | 0.008 mg/plate | TA 1538 | 0.27 | * |
| 0.2 mg/plate | TA 1535 | * | * | 0.0016 mg/pl. | TA 1538 | 0.27 | * |
| 0.04 mg/plate | TA 1535 | * | 0.19 | 0.00032 mg/pl. | TA 1538 | 0.32 | 0.07 |

(act): 5-9 fraction was added

* : calculated value resulted in a negative MUTAR or zero MUTAR

TABLE 6C
MUTAGENIC ACTIVITY RATIO

Substance Assayed: Code #82 Dissolved in: DMSO
Study Number: 81013 Date: 15 July 1981 By: Sauers

| Concentration | Strain | MUTAR (act) | MUTAR | Concentration | Strain | MUTAR (act) | MUTAR |
|----------------|---------|----------------|-------|----------------|---------|----------------|-------|
| 1.0 mg/plate | TA 98 | 0.12 | * | 0.008 mg/plate | TA 1535 | * | 0.06 |
| 0.2 mg/plate | TA 98 | 0.44 | * | 0.0016 mg/pl. | TA 1535 | * | * |
| 0.04 mg/plate | TA 98 | * | * | 0.00032 mg/pl. | TA 1535 | * | 0.13 |
| 0.008 mg/plate | TA 98 | 0.08 | * | | | | |
| 0.0016 mg/pl. | TA 98 | 0.24 | * | 1.0 mg/plate | TA 1537 | * | * |
| 0.00032 mg/pl. | TA 98 | 0.32 | * | 0.2 mg/plate | TA 1537 | 0.31 | * |
| | | | | 0.04 mg/plate | TA 1537 | * | * |
| 1.0 mg/plate | TA 100 | 0.07 | * | 0.008 mg/plate | TA 1537 | 0.31 | * |
| 0.2 mg/plate | TA 100 | 0.18 | * | 0.0016 mg/pl. | TA 1537 | 0.15 | * |
| 0.04 mg/plate | TA 100 | * | * | 0.00032 mg/pl. | TA 1537 | * | * |
| 0.008 mg/plate | TA 100 | 0.23 | * | | | | |
| 0.0016 mg/pl. | TA 100 | * | * | 1.0 mg/plate | TA 1538 | 0.21 | 0.14 |
| 0.00032 mg/pl. | TA 100 | * | 0.07 | 0.2 mg/plate | TA 1538 | 0.21 | 0.07 |
| | | | | 0.04 mg/plate | TA 1538 | 0.27 | * |
| 1.0 mg/plate | TA 1535 | * | * | 0.008 mg/plate | TA 1538 | 0.16 | * |
| 0.2 mg/plate | TA 1535 | * | * | 0.0016 mg/pl. | TA 1538 | 0.16 | * |
| 0.04 mg/plate | TA 1535 | * | 0.26 | 0.00032 mg/pl. | TA 1538 | 0.16 | * |

(act): S-9 fraction was added

* : calculated value resulted in a negative MUTAR or zero MUTAR

TABLE 6D

MUTAGENIC ACTIVITY RATIOSubstance Assayed: Code #87 Dissolved in: DMSOStudy Number: 81013 Date: 15 July 1981 By: Sauers

| Concentration | Strain | MUTAR (act) | MUTAR | Concentration | Strain | MUTAR (act) | MUTAR |
|----------------|---------|----------------|-------|----------------|---------|----------------|-------|
| 1.0 mg/plate | TA 98 | 0.36 | * | 0.008 mg/plate | TA 1535 | 0.09 | 0.13 |
| 0.2 mg/plate | TA 98 | 0.16 | * | 0.0016 mg/pl. | TA 1535 | * | * |
| 0.04 mg/plate | TA 98 | 0.16 | * | 0.00032 mg/pl. | TA 1535 | * | * |
| 0.008 mg/plate | TA 98 | 0.32 | * | | | | |
| 0.0016 mg/pl. | TA 98 | 0.12 | * | 1.0 mg/plate | TA 1537 | 0.31 | * |
| 0.00032 mg/pl. | TA 98 | * | * | 0.2 mg/plate | TA 1537 | 0.46 | 0.15 |
| | | | | 0.04 mg/plate | TA 1537 | 0.31 | * |
| 1.0 mg/plate | TA 100 | 0.17 | * | 0.008 mg/plate | TA 1537 | 0.77 | * |
| 0.2 mg/plate | TA 100 | 0.11 | * | 0.0016 mg/pl. | TA 1537 | * | * |
| 0.04 mg/plate | TA 100 | 0.06 | * | 0.00032 mg/pl. | TA 1537 | 0.31 | * |
| 0.008 mg/plate | TA 100 | 0.04 | 0.05 | | | | |
| 0.0016 mg/pl. | TA 100 | 0.12 | * | 1.0 mg/plate | TA 1538 | 0.05 | * |
| 0.00032 mg/pl. | TA 100 | * | * | 0.2 mg/plate | TA 1538 | 0.37 | * |
| | | | | 0.04 mg/plate | TA 1538 | 0.21 | * |
| 1.0 mg/plate | TA 1535 | * | * | 0.008 mg/plate | TA 1538 | 0.32 | * |
| 0.2 mg/plate | TA 1535 | * | 0.06 | 0.0016 mg/pl. | TA 1538 | 0.11 | * |
| 0.04 mg/plate | TA 1535 | * | 0.39 | 0.00032 mg/pl. | TA 1538 | * | * |

(act): S-9 fraction was added

* : calculated value resulted in a negative MUTAR or zero MUTAR

OFFICIAL DISTRIBUTION LIST

Comander
US Army Medical Research and Development Command
ATTN: SGRD-SI/ Mrs. Madigan
Fort Detrick, Frederick MD 21701

Defense Technical Information Center
ATTN: DTIC-DDA (12 copies)
Cameron Station
Alexandria VA 22314

Director of Defense Research and Engineering
ATTN: Assistant Director, Environmental and
Life Sciences
Washington DC 20301

The Surgeon General
ATTN: DASG-TLO
Washington DC 20314

HQ DA (DASG-ZXA)
WASH DC 20310

Superintendent
Academy of Health Sciences
ATTN: AHS-COM
Fort Sam Houston TX 78234

Assistant Dean
Institute and Research Support
Uniformed Services University of Health Sciences
6917 Arlington Road
Bethesda MD 20014

Commander
US Army Environmental Hygiene Agency
Aberdeen Proving Ground MD 21070

US Army Research Office
ATTN: Chemical and Biological Sciences Division
P.O. Box 1221
Research Triangle Park NC 27709
Biological Sciences Division
Office of Naval Research
Arlington VA 22217

Director of Life Sciences
USAF Office of Scientific Research (AFSC)
Bolling AFB
Washington DC 20332

Director
Walter Reed Army Institute of Research
Washington DC 20012

Commander
US Army Medical Research Institute of Infectious
Diseases
Fort Detrick, Frederick MD 21701

Commander
US Army Research Institute of Environmental
Medicine
Natick MA 01760

Commander
US Army Institute of Surgical Research
Brooke Army Medical Center
Fort Sam Houston TX 78234

Commander
US Army Institute of Dental Research
Washington DC 20012

Commander
US Army Medical Bioengineering
Research and Development Laboratory
Fort Detrick, Frederick MD 21701

Commander
US Army Aeromedical Research Laboratory
Fort Rucker AL 36362

Commander
US Army Biomedical Laboratory
Aberdeen Proving Ground
Edgewood Arsenal MD 21010

Commander
Naval Medical Research Institute
National Naval Medical Center
Bethesda MD 20014

Commander
USAF School of Aerospace Medicine
Aerospace Medical Division
Brooks Air Force Base TX 78235